


ISOTHERM RADIANT WARMER

The Isotherm Radiant Warmer is used to provide safe and efficient warmth to newborns whilst remaining easy and intuitive for clinicians to use.



The Isotherm Radiant Warmer helps to create a warm environment and maintain comfort. With adaptable installation, user-friendly heat control and a clinician-focused design, it provides a reliable radiant warming solution across perinatal care settings.

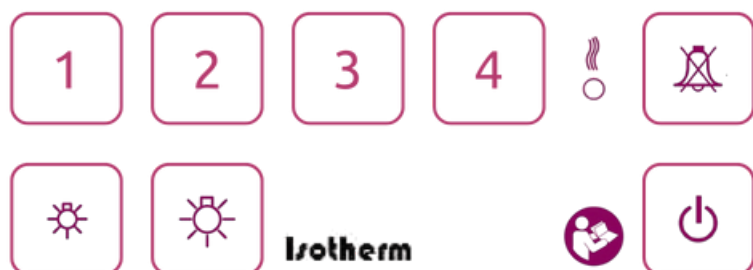

Four clearly defined heat settings make temperature adjustment simple and intuitive, reducing confusion and ensuring consistent care across staff. The factory default settings are as follows:

Level 1: 20% Output | Level 2: 50% Output
Level 3: 75% Output | Level 4: 98% Output



SIMPLE OPERATION

Developed with clinicians in mind, Isotherm's integrated examination light offers both dim and bright settings for clear visibility without disturbing the infant, while the user-friendly interface allows reactive adjustments in busy clinical settings.


Meanwhile, ergonomic side handles let staff adjust the lamp angle safely, avoiding contact with heated surfaces.



ENERGY EFFICIENT

The Isotherm Radiant Warmer delivers exceptional efficiency. Using low-energy LED lighting, it achieves over 95% efficiency for maximum performance with minimal energy use.

By choosing Isotherm, trusts can ensure reliable patient warming while reducing energy use, supporting long-term savings and aligning with the NHS's sustainability commitments.

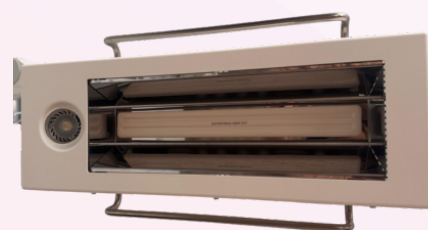

SAFE, CONTROLLED HEAT

No visible infrared light, allowing natural skin tone assessment.

Audio-visual alarms provide clear safety prompts, including a 15-minute review alert to check heat settings.

Durable ceramic heating element ensures uniform, efficient radiant warmth for consistent temperature control in any setting.



VERSATILE MOUNTING OPTIONS

Fixed Wall-mount – ensures consistent warmth is always available at the point of care while keeping equipment securely in place and floors clear for staff movement.

Moving Wall-mount – Provides warmth exactly where it's needed, offering flexible positioning without taking up valuable floor space.

Mobile Stand – allows you to deliver targeted warming support wherever it's required - with a robust adjustable arm and height flexibility.

EVERYTHING YOU NEED TO KNOW

Click the links below to access the Croyde Knowledge Hub

[How to Buy](#)

[Dimensions & Specifications](#)

[Delivery & Implementation](#)

[Warranty & Servicing](#)

[Cleaning Instructions](#)



CROYDE MEDICAL

CONTACT US

Tel: 01367 241194

Email: info@croydemedical.co.uk

[**www.croydemedical.co.uk**](http://www.croydemedical.co.uk)

Croyde Medical, Unit A4 GBP1, Grove Business Park,
Downsview Road, Wantage, OX12 9FA