

SAFEWARM PATIENT WARMING SYSTEM

The Safewarm Patient Warming System provides regulated heat to help patients maintain a comfortable body temperature in operating rooms, lowering their risk of hypothermia.



The Safewarm Mattress & Blanket can reach 37°C in just 5-7 minutes, offering fast, homogenous and reliable warming to patients that need their body temperature regulated.

There is an accuracy of 0.1°C when controlling the mattress and blanket, offering additional sensitivity and flexibility to health professionals when determining the best temperature for their patient.

Built-in sensors mean the mattress or blanket will automatically cut-off heat if it reaches a temperature of 42°C, ensuring it remains safe by avoiding overheating.



INFECTION CONTROL

By using conductive heat rather than forced air, they avoid circulating air around the patient, helping to protect sterile fields and reduce the risk of airborne contamination.


This makes Safewarm an excellent choice for orthopaedic and plastic or reconstructive procedures.



MULTI-CHANNEL CONTROL

The innovative control unit comes with three channels, meaning up to three mattresses and/or blankets can be operated at the same time.

This offers enhanced flexibility over patient warming configurations.

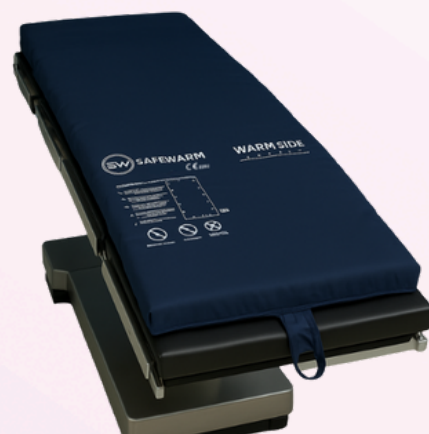

MINIMAL UPKEEP

The mattresses and blankets can be easily cleaned using hospital-approved disinfectants.

There are also reusable, washable covers available to further enhance the ease of cleaning and speed up turnover between patients.



DESIGNED TO BE SAFE

Independent microprocessors within the mattress, blanket and controller will set off an audio-visual alarm if any faults are detected, allowing for prompt intervention.

The inner mattress is made from a carbon fibre material that allows X-rays to pass through, enabling patient temperatures to be maintained during imaging procedures.

Heat is generated with a protective low voltage of 24V, mitigating any risk of electric shock.

EVERYTHING YOU NEED TO KNOW

Click the links below to access the Croyde Knowledge Hub

[How to Buy](#)

[Dimensions & Specifications](#)

[Delivery & Implementation](#)

[Warranty & Servicing](#)

[Cleaning Instructions](#)



CROYDE MEDICAL

Maternity Equipment Specialists

CONTACT US

Tel: 01367 241194

Email: info@croydemedical.co.uk

[**www.croydemedical.co.uk**](http://www.croydemedical.co.uk)

Croyde Medical, Unit 7, Regal Way, Faringdon, Oxon, SN7 7BX